

Aluminum Lightweight Wheeled Stretcher

Lightweight Wheeled Stretcher (High-Low Model)

Crank
Handle
Height
Adjustable

Fixed
Stretcher


Easy to operate for all users.



AL-UD-1 (4-Wheel Double Lock Model)

(Tax Included) **¥330,088**

- Material/Aluminum
- Dimension/Total Width 590 x Total Length 1,900 x Total Height 525~910mm
- Casters/φ150mm Double Lock
- Weight/40.2kg
- Method of Height Adjustment/Crank Handle
- Standard Features/Mattress, Irrigator Stand
- Mattress Dimensions/Total Width 500 x Total Length 1,690 x Total Height 50mm

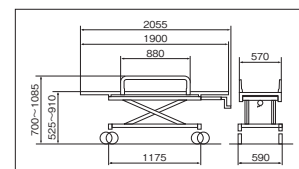


Crank
Handle
Height
Adjustable

Fixed
Stretcher

AL-TUD-1 (Total Lock Model)

(Tax Included) **¥431,244**

- Material/Aluminum
- Dimension/Total Width 590 x Total Length 1,900 x Total Height 525~910mm
- Casters/φ150mm Total Lock U-Shaped Model
- Weight/42.8kg
- Method of Height Adjustment/Crank Handle
- Standard Features/Mattress, Irrigator Stand
- Mattress Dimensions/Total Width 500 x Total Length 1,690 x Total Height 50mm

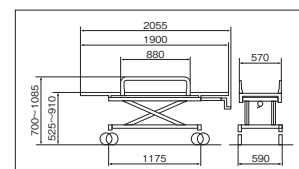

Crank
Handle
Height
Adjustable


Fixed
Stretcher

* Pediatric Model of **AL-UD-1**

AL-UD-1S (4-Wheel Double Lock Model - Pediatric)


(Tax Included) **¥408,617**

- Material/Aluminum
- Dimension/Total Width 500 x Total Length 1,500 x Total Height 610~920mm
- Casters/φ150mm Double Lock
- Weight/34.0kg
- Method of Height Adjustment/Crank Handle
- Standard Features/Mattress, Irrigator Stand
- Mattress Dimensions/Total Width 390 x Total Length 1,300 x Total Height 50mm

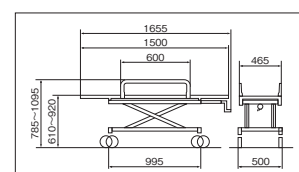

What are double lock casters • total lock casters?

Double Lock Casters

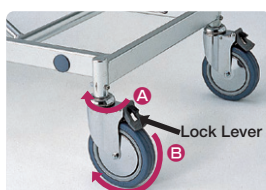

Standard
φ150mm

Lowering the lock lever simultaneously locks swivel movement **A** and rotational movement **B** of the casters.

Total Lock U-Shaped Model



Standard
φ150mm

The following 2 locking patterns are engaged by raising and lowering the lock lever.

① Simultaneously locks movements **A** and **B** for all 4 wheels.

② Locks movement **A** of only 1 caster.

This stretcher is useful for transporting patients in straight lines over long distances, e.g. in hospital hallways.

Standard Feature

Irrigator Stand Set

